


Face Shields

How to Remove the Liner

Our face shields come with a protective plastic film as a liner to prevent scratches to the shields during packing and shipping. Remove the liner before using the product.

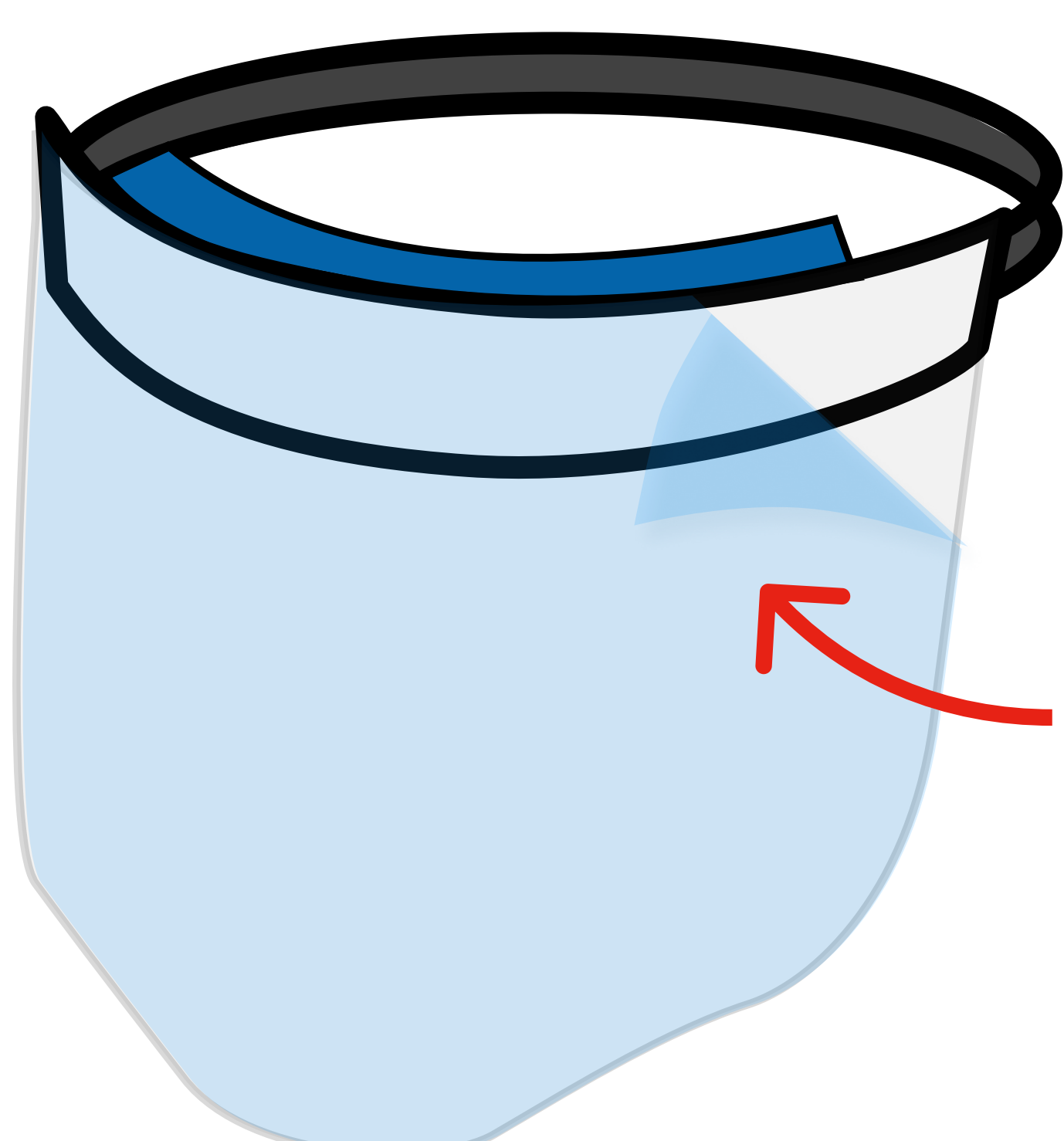

- 1 With a piece of tape (Scotch or another tape with a good grip), attach part of the tape to the corner of the plastic film.
- 2 Hold the non-attached part of the tape and pull backwards to peel the film away from the shield.
- 3 If you are having difficulty, try holding down the shield area just below the tape while you peel. You can also use tweezers to remove the film from its corners.

How to Wear the Face Shield



- 1 Remove the liner film before wearing the face shield.
- 2 Stretch the elastic around your head so that the foam sits on your forehead.
- 3 Adjust as necessary.

Please note that this product is suitable for individual use only and should not be shared with anybody else.

How to Clean the Face Shield

It is important to clean the face shield after each use. We recommend cleaning and sanitizing the face shield in the following ways.

- 1 Wash both sides of the face shield with soap and warm water.
- 2 Apply sanitizer on a paper towel and use it to wipe down the face shield. Do not use aggressive cleaning agents.
- 3 Store the face shield in a safe, clean, and dry location. You can also store it in a large Ziploc bag.
- 4 Dispose of the face shield after seven days of use.